

SPOKE GUARDS

- p) Earth
- q) Lion
- r) Flamingo



- s) Orange
- t) Lightning
- u) Girl



Spoke guards with custom design What supplies Mobilex?

Mobilex provides personalized spoke guards in sizes 20", 22" and 24". Other dimensions depends on quantities and can be offered on request.

Spoke guards are silk screen printed (serigraph).

Spoke guards are made of transparent PETG, are slightly curved and have cutouts for the 4 or 6 fastening segments of the push rims. They are fastened with 6 clips on the spokes of the wheelchair wheel. The spoke guards have the required holes for the fastening clips, and the required clips are included. Shape and material correspond to the printed spoke guards from the Mobilex standard program.

Mobilex takes care of the whole process of the production, the printing and the delivery.

Spoke guards with custom design What will we need from you?

Formats:

Under ideal circumstances you provide us with the desired drawing or logo already processed as vector graphics (.svg, .eps, .ai) because this enables us to make the layout easily.

If your graphical material is not available as vector graphics, we also accept .tif, .png or .jpg formats with a resolution of 300 DPI at a size of approximately: 200 x 200 mm.

Colours:

The colours should be defined with Pantone or HKS colour system codes. A printing film is created for each colour, then, colours are placed separately one after the other on the spokes guards. Because of this we need for each colour the exact indication of the Pantone or HKS colour codes.




SPOKE GUARDS



Händler

SPOKE GUARDS
SPOKE GUARDS

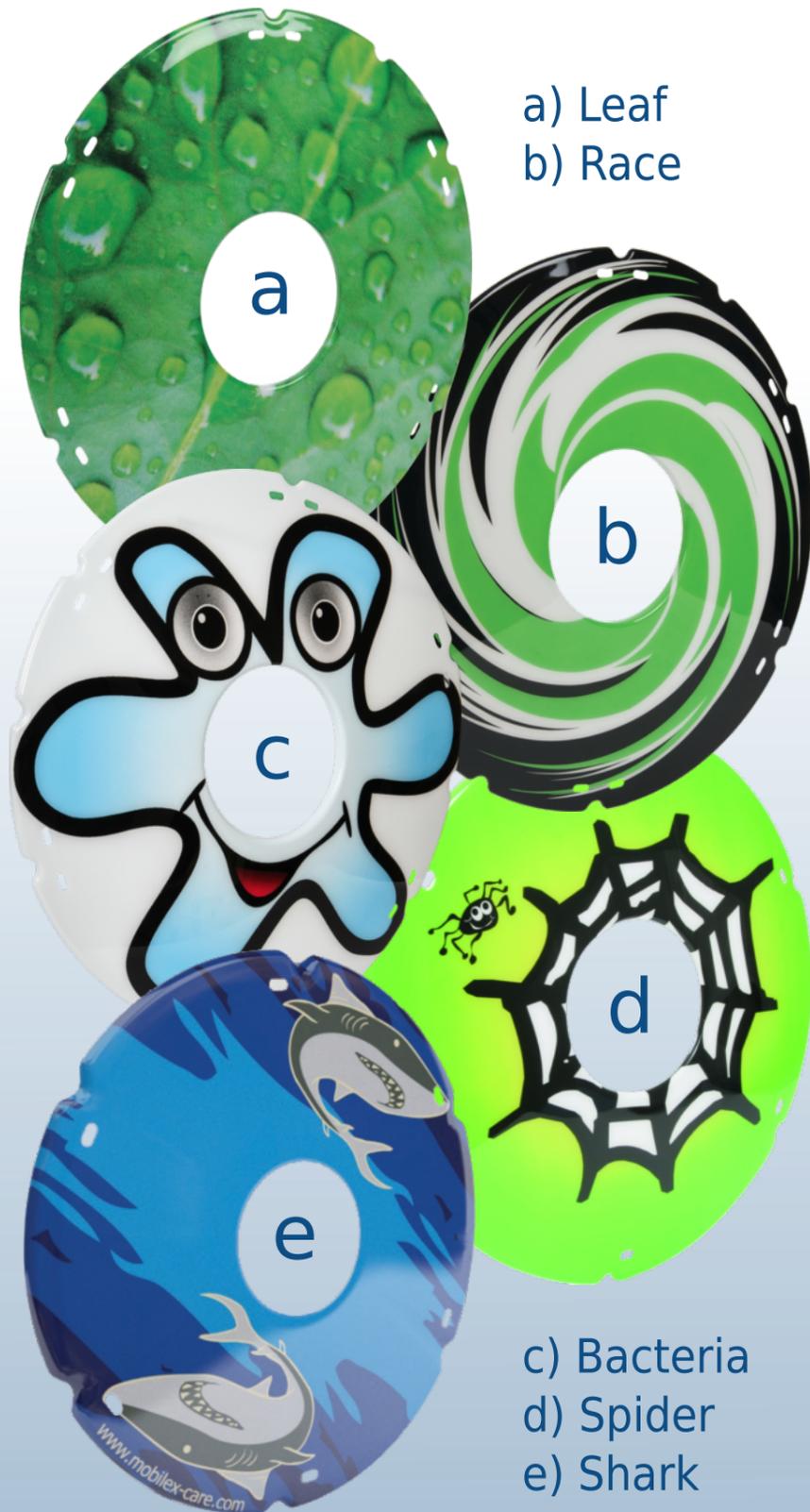

SPOKE GUARDS
SPOKE GUARDS



SPOKE GUARDS
SPOKE GUARDS

