


TIGER

OUTDOOR ROLLATOR



Mobilex
a better life in motion

Rollator Tiger * with rotatable armrests	Standard	Outdoor
	Item Nr. 312201 (grey) Item Nr. 312211 (red) Item Nr. 312202 (brown) *	Item Nr. 312205 (grey) * Item Nr. 312206 (brown) *
Height of the forearm support (adjustable)	90 cm - 115 cm	87 - 110 cm
Recommended user heights	150 - 200 cm	140 - 185 cm
Seat height	62 cm	62 cm
Total width opened / folded	62 cm / 22,5 cm	71 cm / 31 cm
Seat width / Inner width between push handles	46 cm / 36 cm	56 cm / 46 cm
Length (depth)	66 cm	77 cm
Weight (+ bag 380 g)	9,5 kg	12 kg
Max. load (rollator / bag)	150 kg / 5 kg	200 kg / 5 kg
Brakes	2 lockable with parking function	
Wheel sizes	200 x 35 mm	285/235 x 44 mm

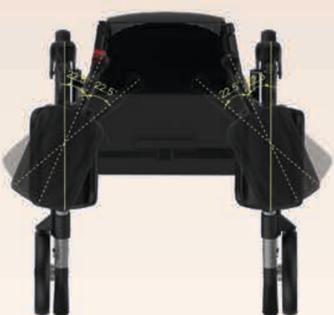

The Mobilex Arthritis Rollator Tiger Outdoor is based on the elegant Lion rollator and therefore has all its features and functions. As an Arthritis Rollator it has soft and comfortable forearm supports. These can be adjusted in depth, angle and height and can be turned inwards in 2 angles by 22.5 degrees each. The rollator is extremely stable and easy to manoeuvre. It has large, smooth-running wheels and has smooth-running and lockable brakes. In the driving position it is secured against unintentional unlocking. As standard, the rollator has a transport lock and a removable shopping bag. A wide range of accessories are available. The Tiger Outdoor model is approved for outdoor use with a user weight of 200 kg.

