

TRANSFER PRODUCTS

MOBILEX SELECTION OF TRANSFER PRODUCTS



TRANSFER BOARD FROM PVC

Item nr. 240505

The Mobiglide transfer board is made of 6 mm thick white PVC. The low frictional resistance of the Mobiglide transfer boards makes transferring to or from the wheelchair much easier. Anti-slip strips are available as accessories, which can be attached to the underside of the board to prevent it from sliding away. The transfer board is loadable up to 150 kg with a maximum intermediate distance of 20 cm.

The PVC board is available in the size 800 × 245 mm.

TRANSFER BOARD FROM POLYPROPYLENE

Item nr. 240501, 240502 & 240507

The Mobiglide transfer boards are made of 6 mm thick white translucent polypropylene. The low frictional resistance of the Mobiglide transfer boards makes transferring to or from the wheelchair much easier. Anti-slip strips are available as accessories, which can be attached to the underside of the board to prevent it from sliding away. The transfer board is loadable up to 150 kg with a maximum intermediate distance of 20 cm.


The polypropylene boards are available in three sizes: 650 × 245 mm, 800 × 245 mm and a curved form 800 × 280 mm.



MOBIGLIDE PATIENT TRANSFER BOARD

Item nr. 240701

The Mobiglide patient transfer board made of break-proof 6 mm thick HDPE has a high sliding ability on the upper side and helps when transferring patients from bed to bed. In hospitals or nursing homes it can also be used in oncology for X-rays as the material allows X-rays to pass through. Anti-slip strips are available as an accessory that can be attached to the underside of the board to prevent it from sliding away. The Mobiglide patient transfer board has a size of 1500 × 700 mm. The load capacity is 200 kg with a max. distance of 20 cm and a max. height difference of 10 cm.

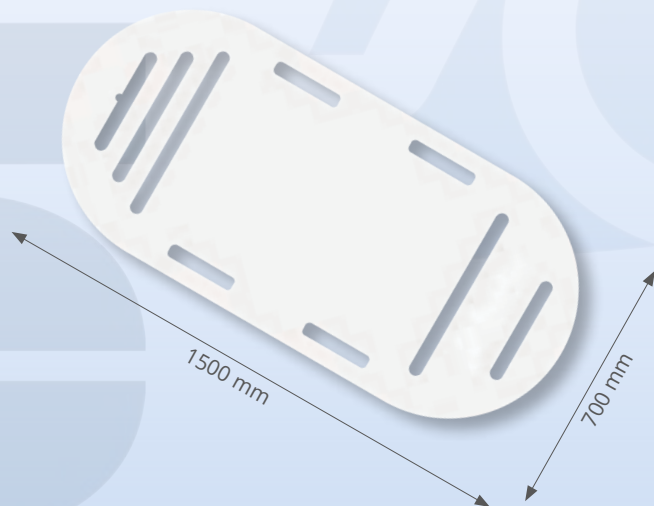

SUPPORT BELT

Item nr. 278046 & 278047

The Mobiglide support belt helps caregivers to lift patients out of bed, chair or wheelchair and to travel short distances with them. This can be a supportive aid when getting up from sitting or vice versa, e.g. from the wheelchair or bed, the car, the toilet, etc. Furthermore, the support belt can be used to support the patient during gait training.

Sizes: 117 × 14 cm (M) or 140 × 14 cm (L)
Load capacity: 180 kg





FLAT SLIDING SHEET

Item nr. 278015 & 278016

The Mobiglide flat sliding sheets reduce the frictional forces when transferring or moving people in bed. The sliding sheet is to be used in pairs.

Size: 150 × 100 cm or 200 × 100 cm
Load capacity: 250 kg

TUBE SLIDE SHEET

Item nr. 278011, 278012 & 278013

The Mobiglide tube slide sheets offer extremely low frictional resistance on the inside and outside, thus minimizing friction when transferring patients. The tube slide sheets are available in 3 sizes.

Size: 70 x 68 cm, 140 x 68 cm and 200 x 68 cm
Load capacity: 250 kg



SLIDE SHEET WITH HANDLES

Item nr. 278023

The Mobiglide slide sheet with handles makes positioning of sitting or lying patients easier than well as transferring a lying patient from a bed to a stretcher. The use of two tube slide sheets makes it easier to place the lifting sling under the patient.

Size: 100 × 200 cm
Load capacity: 250 kg

ONE-WAY SLIDE MAT

Item nr. 278031

The Mobiglide one-way slide mat prevents the user from gliding forward in the wheelchair. The internal velour can only slide in one direction and prevents the user from slipping unintentionally forward. The sliding mat also facilitates the user to transfer back to his original seating position.

Size: 45 × 50 cm
Load capacity: 250 kg



SLIDING GLOVES

Item nr. 278041

The Mobiglide sliding gloves make it easier to slide under the patient, e.g. to check the filling level of an anti-decubitus product. They also make it easier to carry out positioning or repositioning.

Size: 50 × 20 cm



STANDARD SWIVEL SEAT CUSHION

Item nr. 278051 & 278052

The standard swivel seat cushion consists of two textile parts joined together. The smooth inner sides make it easier to turn. The upper part consists of a dark blue soft surface, the lower part has a non-slip underside. The cushion is intended to be used as a seat pad in a car, armchair or chair to facilitate turning and transferring of the user.

Diameter: 45 or 50 cm
Max. Load capacity: 100 kg



MAGIC-TURN SWIVEL SEAT CUSHION

Item nr. 278055 & 278056

The Magic Turn swivel seat cushion consists of two textile parts joined together. The "Magic Turn" model offers greater comfort due to the thicker padding. The insides are slide-coated for excellent turning properties. The bottom part has a non-slip underside. The cushion is designed to be used as a seat base in a car, armchair or chair to facilitate turning and transferring of the user.

Diameter: 37 or 44 cm
Max. Load capacity: 130 kg



MAGIC-TURN SWIVEL STANDING CUSHION

Item nr. 278058

The Magic-Turn swivel standing cushion consists of two textile parts connected to each other. The insides are slide-coated to make it easier to turn. The upper part has a soft surface, the lower part is equipped with a non-slip underside. The swivel cushion is intended as a transfer aid, e.g. from bed to chair, etc., and supports turning on its own axis. Can also be used in wet areas as both the top and bottom have an anti-slip coating.

Diameter: 37 cm
Max. Load capacity: 130 kg

